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DEVICE FAMILY: Contact lens Hydrogel reusable https://www.lens55.com/descargas.html Do not flush in toilet **BC** Base Curve **DIA** Diameter **Pwr** Power

The expiry date indicated on the contact lens packaging is the date on which you can still open the blister and wear the contact lenses for 30 days.

Models #	**Material** 18369-1:2017	**Replacement**
L551.1 L551.2 L551.3 L551.4 L551.5 L551.6	Ocufilcon D 4(19)[55%] \| Methafilcon A 4(19)[55%]	Monthly/Daily Wear<24 hours Discard after 30 days open
Sil2.1 Sil2.2 Sil2.3 Sil2.4	Filcon 5 B (50)[58%]	Monthly/Daily Wear<24 hours Discard after 30 days open
PC3.1 PC3.2 PC3.3 PC3.4	Filcon 5 B (50)[60%]	Monthly/Daily Wear/Flexible wear Discard after 30 days open
Col4.1 Col4.2 Col4.3 Col4.4	Deltafilcon A 1 (10)[43%]	Monthly/Daily Wear<24 hours Discard after 30 days open

EN - INSTRUCTIONS FOR USE FOR CONTACT LENS HANDLING, INSERTION AND REMOVAL

INFORMATION LEAFLET, TO BE PROVIDED TO WEARERS OF CONTACT LENSES BY THEIR EYE CARE PROFFESIONAL (optometrist/ophthalmologist)

THE PURPOSE OF THIS INSTRUCTIONS (MANUAL) IS TO GUIDE THE WEARER HOW TO PROPERLY WEAR AND CARE FOR CONTACT LENSES. FOLLOW THE INSTRUCTIONS OF YOUR EYE CARE PROFESSIONAL, BEFORE THE USE OF CONTACT LENSES, PLEASE READ CAREFULLY AND KEEP THIS INFORMATION FOR FUTURE USE.

Contact lenses should always be fitted by an eye care professional. Servilens manufactures different contact lenses and, therefore, the eye care professional is qualified to recommend the contact lens that best suits your needs in terms of parameters, materials, replacement and best use for your vision.
The traceability of contact lenses is guaranteed by an eye care professional.

INTENDED PURPOSE / INTENDED USE

Servilens contact lenses are designed to be worn on the front of the eye and are mainly intended to correct defects on visual function.

Indications

Servilens contact lenses are a medical device indicated for correcting eye disorders (myopia [short-sightedness], hyperopia [long-sightedness], presbyopia and astigmatism) in healthy eyes for daily wear, as long as it is comfortable.
Do not wear the lenses while sleeping, unless advised by an eye care professional.

INTENDED USERS/TARGET USER GROUP

People with healthy eyes who need a correction of eye disorders. There is no age restriction. The appropriate age should be determined by an eye care professional.
People with corneal dystrophies, corneal ulcers, inflammation of cornea or the eye, chronic dacryocystitis, chronic blepharitis, chronic conjunctivitis or allergic eye conditions, must not use contact lenses.

INSTRUCTIONS FOR HANDLING YOUR CONTACT LENSES

See in figures 1 – 5.
Your eye care professional should provide you with detailed contact lenses insertion and removal instructions, as well as information on the safe use of your lenses.
Wash your hands with soap and water; this decreases the chance of an eye infection.
Dry your hands with tissues or a lint-free cloth so that particles don't wind up floating around in one or both of your eyes.
Always wear the right contact lens on the right eye first and then wear the left one. Always remove the contact lens from the right eye first and once the right lens has been stored and cleaned in the case, remove the left lens. This keeps you from mixing up the right-eye and left-eye contacts, which typically have different power or measurements.
Do not use sharp objects or your fingernails to handle your contact lenses.

BLISTER (PP material + Aluminum) - CONTACT LENS SOLUTION

The soft contact lens is immersed in buffered 0.9% saline solution (sodium chloride, boracic acid, borax, purified water). Gently shake the blister pack before opening. Peel back the foil lid. Pour the lens into the palm of your hand or if necessary, gently slide the lens out of the container, or open as directed by your eye care professional.

CONTACT LENSES MAINTENANCE

Use soft contact lens care solution (All in One or Peroxide)
After removing contact lenses, always clean and disinfect them according to the instructions on the care solution.

IT IS ESSENTIAL THE CORRECT CARE OF YOUR CONTACT LENSES. Your eye care professional should choose the lens care system most suitable for you. Do not use heat disinfection. NEVER allow contact lenses to come in contact with non-sterile water (including tap water) as microbial contamination can occur, which may lead to permanent eye damage. Do not change the contact lens maintenance without contacting with your eye care professional. Ask your eye care professional about wearing contact lenses during sporting activities. Remove lenses when exposed to noxious or irritating vapors.
If the contact lens is kept and not worn over one week, re-cleaning the lens is necessary before wearing the contact lens and contact with your eye care professional.
Contact lenses should be stored at room temperature, not in the fridge. Avoid extreme temperatures.

CONTACT LENS CASE

Lens case is the important source of microbial pollution. In order to prevent the eyes infection, the lens case needs to be cleansed, washed after everyday use. Please regularly change the lens case according to the time recommended by your eye care professional.

CONTACT LENSES CONTRAINDICATIONS

Do not use contact lenses in corneal dystrophies, corneal ulcers, inflammation of cornea or the eye, chronic dacryocystitis, chronic blepharitis, chronic conjunctivitis, allergic eye conditions.
Do not use contact lenses if cannot inspect your eyes regularly.

CONTACT LENS SIDE-EFFECTS

Contact lenses can cause permanent damage if used improperly. When you are wearing contact lenses, check your eyes every day. Your eyes should look well, you should feel comfortable with them, and your vision should be clear. If it is not, remove your lenses immediately and contact your eye care professional.
Unusual symptoms, if ignored, can lead to serious complications.
In the event of a serious incident involving the product, the eye care professional shall be informed and shall be responsible for communicating the problem to the manufacturer and the competent health authorities.

CONTACT LENSES - IMPORTANT PRECAUTIONS

If there are any problems with the lens (dry out, deformation, discoloration, foreign bodies, damage, scratches or dirt), discard the contact lens and consult your care professional.
Your eye care professional will recommend a lens wearing, cleaning and replacement regime.
If your eyes become red or irritated or you experience any pain or unexpected change in vision, remove your lenses immediately and consult your eye care professional.
Avoid contaminating contact lenses with make-up, aerosol sprays or wearing them near noxious fumes.
Replace damaged, contaminated or scratched lenses immediately.
Avoid long nails. A scratched eye is vulnerable to infection.
Keep your contact lenses far from the place that child can touch.
Only use your contact lenses for water sports if watertight goggles are worn.
Consult your eye care professional about using lenses while performing activities that could be considered dangerous.
Contact with unsterile water (e.g., from the tap, bottled, bath, pool, or water sports) increases the risk of microbial infection.
If you wear contact lenses frequently and satisfactorily and within the first few minutes of wearing your contact lenses feel irritant symptoms such as blurred vision, redness of the eyes, photophobia, tearing, etc., you should immediately remove the lens and look for the cause. Reasons that may cause the above-mentioned problems are that the lens is not positioned in the middle of cornea; the lens is contaminated and/or damaged; the right and left lens are out of order or other reasons. If the symptoms continue, you should consult your eye care professional.
If the symptoms described above occur after several hours of comfortable contact lens wear, you should immediately remove the lens and consult your eye care professional.


CONTACT LENSES WARNINGS

Drive your motor vehicle after completely adapting to wearing contact lenses. Consult your eye care professional.
Do not share your contact lenses with anyone else.
Do not wear your contact lenses during illness, such as cold, fever, etc.
Do not wear your lenses when using eye drops (unless advised to do so by your eye care professional).
Do not use saliva to wet your lenses.

Do not wear your contact lenses beyond the period recommended by your eye care professional.
Do not use contact lenses if the package is opened or damaged.
Do not rinse your contact lenses case with tap water.
Do not use the contact lenses if the lenses are damaged.
Do not sleep with daily contact lenses.

WEARING SCHEDULE
Follow the wearing schedule provided by your eye care professional.
Wear the contact lenses during the time indicated by your eye care professional. Check your vision periodically with your eye care professional and follow the instructions use.
The lenses are suitable for daily wear. Do not wear the lenses while sleeping, unless your eye care professional has advised it is safe to do so for up to 6 days. Clinical evidence has shown that if contact lenses are worn while sleeping, the risk of serious adverse effects increases considerably. Strictly follow the recommendations of your eye care professional.

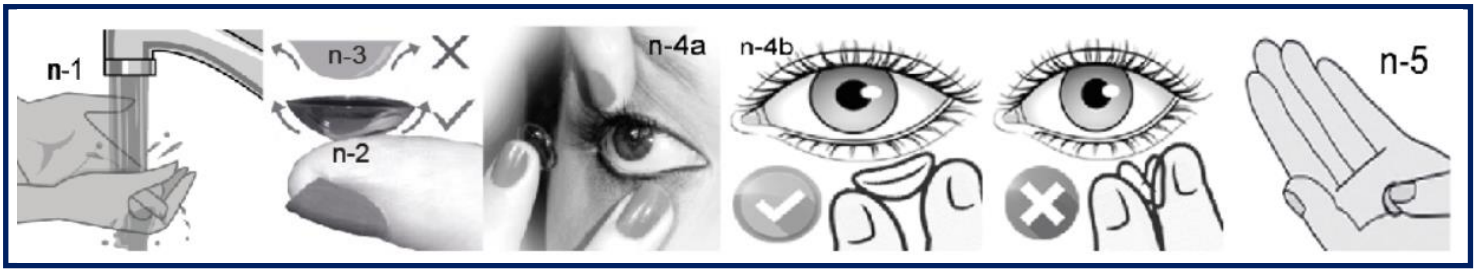Sterility: soft contact lenses in blister are sterilized by steam (sterility is only guaranteed in the closed and undamaged sealed package)
No special conditions are required for the disposal of the blister and contact lens. Do not flush in toilets.
Product name/parameters: see packing label. **Material name**: see packing label **Administration**: suitable for eye use.

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